

Premisut Catgut Plain

Premisut

Description

Premisut Catgut Plain is a sterile absorbable surgical suture composed of purified connective tissue (collagen) derived from the sub-mucous layer of sheep (bovine) intestine. Premisut Catgut Plain is available as Plain and in a range of gauge sizes and lengths, non-needled or attached to stainless steel needles of varying types and sizes. Premisut Catgut Plain is packed in conditioning fluid, the constituents of which are isopropanol, water, Sodium Benzoate and Ethylene oxide.

Indications

Premisut Catgut Plain is indicated for use in soft tissue approximation and or ligation, including use in plastic & ophthalmic procedures but not for use in cardiovascular and neurological tissues.

Performance


When Premisut Catgut Plain suture is implanted, a moderate tissue inflammation occurs which is characteristic of a foreign body response. Loss of tensile strength and loss of suture mass follows as the proteolytic enzymatic digestive process proteolyse the surgical catgut. Due to the inherent variability of a natural material, figures for strength loss and absorption following implantation are for guidance only. Effective strength is usually lost between 10 to 12 days in case of Premisut Catgut Plain. Absorption is essentially complete by 90-100 days.

Total Absorption

~90~100 days

Color

Ivory

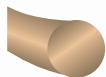

Sizes

5/0 to 2 USP

Usage Areas

General surgery, gastrointestinal, gynecology, obstetrics, urology, plastic, ophthalmic.

Sterilization

Ethylene Oxide

Shelf Life

5 years

Packaging

Box of 36pcs



Premisut CATGUT PLAIN

● ½ Circle Round bodied TAPERPOINT 

NEEDLE LENGTH MM		USP	EP	COLOUR	THREAD LENGTH CM	ATLAS REF	PCS/BOX
HR 16		5/0	1.5	ivory	75	AMCP2HR1650	36
		4/0	2	ivory	75	AMCP2HR1640	36
HR 20		4/0	2	ivory	75	AMCP2HR2040	36
		3/0	3	ivory	75	AMCP2HR2030	36
		2/0	3.5	ivory	75	AMCP2HR2020	36
• RIBBED							
HR 26		4/0	2	ivory	75	AMCP2HR2640	36
		3/0	3	ivory	75	AMCP2HR2630	36
		2/0	3.5	ivory	75	AMCP2HR2620	36
		0	4	ivory	75	AMCP2HR2610	36
• RIBBED							
HR 30		3/0	3	ivory	75	AMCP2HR3030	36
• RIBBED							
HR 37		3/0	3	ivory	75	AMCP2HR3730	36
		2/0	3.5	ivory	75	AMCP2HR3720	36
		0	4	ivory	75	AMCP2HR3710	36
		1	5	ivory	75	AMCP2HR3701	36
• RIBBED							

CATGUT PLAIN

● ½ Circle Round bodied TAPERPOINT 

NEEDLE LENGTH MM		USP	EP	COLOUR	THREAD LENGTH CM	ATLAS REF	PCS/BOX
HR 40		1	5	ivory	90	AMCP3HR4001	36
• RIBBED							
HR 45		0	4	ivory	75	AMCP2HR4510	36
		1	5	ivory	75	AMCP2HR4501	36
		2	6	ivory	75	AMCP2HR4502	36
• RIBBED							
HR 50		2/0	3.5	ivory	75	AMCP2HR5020	36
• RIBBED							

CATGUT PLAIN

▼ 3/8 Circle Reverse Cutting 

NEEDLE LENGTH MM		USP	EP	COLOUR	THREAD LENGTH CM	ATLAS REF	PCS/BOX
DS 12		5/0	1.5	ivory	45	AMCP1DS1250	36
DS 16		6/0	1	ivory	45	AMCP1DS1660	36
		5/0	1.5	ivory	75	AMCP2DS1650	36
		4/0	2	ivory	75	AMCP2DS1640	36
		3/0	3	ivory	75	AMCP2DS1630	36
DS 19		4/0	2	ivory	75	AMCP2DS1940	36
		3/0	3	ivory	75	AMCP2DS1930	36
		2/0	3.5	ivory	75	AMCP2DS1920	36
• RIBBED							
DS 24		4/0	2	ivory	75	AMCP2DS2440	36
		3/0	3	ivory	75	AMCP2DS2430	36
		2/0	3.5	ivory	75	AMCP2DS2420	36
		0	4	ivory	75	AMCP2DS2410	36
• RIBBED							
DS 37		2/0	3.5	ivory	75	AMCP2DS3720	36
		0	4	ivory	75	AMCP2DS3710	36
• RIBBED							

CATGUT PLAIN

Without needle

NEEDLE LENGTH MM		USP	EP	COLOUR	THREAD LENGTH CM	ATLAS REF	PCS/BOX
	Without needle 150cm	4/0	2	ivory	150	AMCP-15040	36
		3/0	3	ivory	150	AMCP-15030	36
		2/0	3.5	ivory	150	AMCP-15020	36
		0	4	ivory	150	AMCP-15010	36
		1	5	ivory	150	AMCP-15001	36
		2	6	ivory	150	AMCP-15002	36